
Mölnlycke®
Tubigrip®

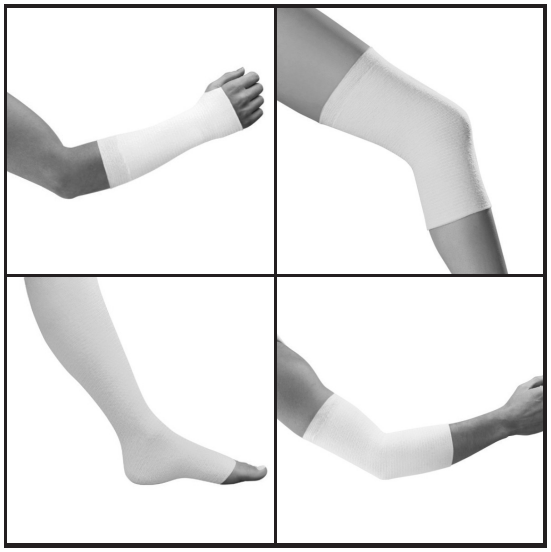en Elasticated Tubular Bandage. Provides light compression. **Precautions:** Do not use Tubigrip if known hypersensitivity to the ingoing materials. If the condition worsens, stop use and consult a health care professional. If any serious incident has occurred in relation to the use of Tubigrip, it should be reported to Mölnlycke Health Care and to your local competent authority. **Instructions for use:** Consult size guide. Cut the bandage to twice the length of the affected limb, if required make 2 slits for a thumb hole. Pull over the limb, then double back. Single patient multiple use. Keep clean. If soiled, wash or replace bandage. Washable up to 5 times, as recommended. Disposal according to local environmental procedures. **Material Content:** Cotton, Elastodiene and Polyester. **Caution: This product contains natural rubber latex which may cause allergic reactions, including anaphylactic responses.**

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LATEX

Caution: This Product Contains Natural Rubber Latex Which May Cause Allergic Reactions, including anaphylactic responses.

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UK Responsible Person:
Mölnlycke Health Care Ltd,
Medlock Street
Oldham
Lancashire
OL1 3HS
United Kingdom

Australian Sponsor Address:
Mölnlycke Health Care Pty. Ltd
12 Narabang Way
Belrose, NSW 2085
Australia

Made in UK

Toll free number :
USA 1-800-882-4582
Canada 1-800-494-5134

www.molnlycke.com
www.molnlycke.com/symbols

CH REP
Mölnlycke Health Care AG
Brandstrasse 24
8952 Schlieren

Manufacturer / Fabricant
Mölnlycke Health Care AB
Entreprenörsstråket 21, SE-431 53 Mölndal, Sweden

The Mölnlycke and Tubigrip trademarks, names and logos are registered globally to one or more of the Mölnlycke Health Care group of companies.

Revised 2026-07
Master PD-774502 rev 3

327200-46
PD-828094 rev. 1