Mölnlycke®
Setopress®


gb P.E.C MEDIUM / HIGH COMPRESSION BANDAGE

Setopress is a polyamide, elastane and cotton compression (medium/high) extensible (“PECCHE”) bandage.


Indications for use
Setopress is for the treatment of venous leg ulcers.

Instructions for use
To ensure correct application, a simple visual guide is printed on both sides of the bandage. The visual guide consists of coloured rectangles which become squares when the correct extension is reached. On one side of the bandage are brown rectangles indicating high compression and on the other side are green rectangles indicating moderate compression.


Before Extension

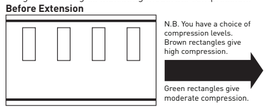Correctly extended - DO NOT OVERSTRETCH


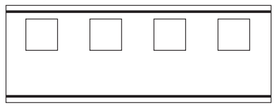
		Applied spirally with 50% overlap			
		Limb circumference (cm/in.)			
	Extension +/- 10%	<18cm	18 – 26cm/ 7 – 10.2in.	27 – 35cm/ 10.6 – 13.8in.	36 – 50cm/ 14.2 – 19.7in.
Approx sub bandage pressure	70% (green square)	Not suitable	27 -38 mmHg/ 1 – 1.5in.Hg	20 – 26 mmHg/ 0.8 – 1.0in.Hg	Not suitable
	100% (brown square)	Not suitable	35 -50 mmHg/ 1.4 – 2in.Hg	27 – 34 mmHg/ 1 – 1.3in.Hg	19 - 26 mmHg/ 0.7 – 1in.Hg

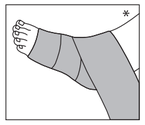
*A padding layer should be used prior to application

Commence bandaging at the base of the toes and enclose heel with NO EXTENSION...

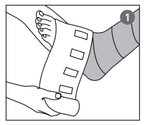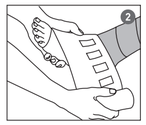
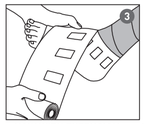
1. Apply around foot starting at base of toes (no extension)
2. Make two turns around the ankle to enclose heel (no extension)
3. Return to base of toes and make one turn around foot (no extension)

...NOW WITH EXTENSION continue up the leg in an ascending spiral

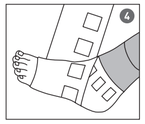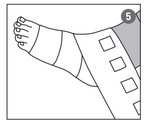
4. Now apply full extension (rectangles becomes squares)
5. Apply with 50% overlap (full extension)
6. Continue to apply as an ascending spiral (full extension)


N.B. It is important to use the same indicator for consistent extension throughout the application.

...finally secure bandage just below the knee

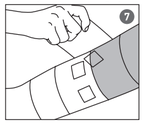7. Stop bandaging just below the knee. Cut off any excess bandage.
8. Secure.

Setopress may be applied over a primary wound dressing.

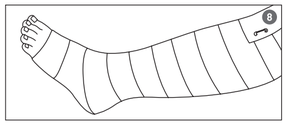Contraindication(s)
Do not use figure of eight bandaging.
Do not apply if arterial insufficiency is present or suspected.
Do not use on limbs with a circumference of less than 18cm/7in.

Precautions
Do not use Setopress on patient and/or user with known hypersensitivity to the ingoing materials/components of the product.
The effective circumference of the ankle when reading off the sizing table must take into account any padding layer.
Setopress should only be applied by a trained bandager.
Do not use bandage clips.
This product must be used under the direction of a healthcare professional.
Bandage should be changed in line with dressing changes.
Discard and replace if soiled with bodily fluids.

Fibre content
Polyamide, Elastane, Cotton.

Special storage conditions and handling conditions
Disposal should be handled according to local environmental procedures.

Other information
If any serious incident has occurred in relation to the use of Setopress, it should be reported to Mölnlycke Health Care and to your local competent authority.
The Mölnlycke and Setopress trademarks, names and logos are registered globally to one or more of the Mölnlycke Health Care Group of companies.
Size 10cm x 3.5m (relaxed roll)

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