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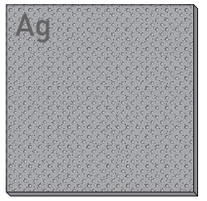
www.molnlycke.com

Revised 2025-02
Master PD-431043 rev. 10
40326-45
Made in Finland
PD-755177 rev 1

en Mepilex® Transfer Ag

Antimicrobial soft silicone exudate transfer dressing

℞ only **Caution: Federal [US] Law restricts this device to sale by or on the order of a physician [or properly licensed practitioner].**

Product description

Mepilex Transfer Ag is a soft silicone wound contact layer that absorbs and transfers exudate, maintains a moist wound environment and has antimicrobial properties.

Mepilex Transfer Ag contains silver sulfate which acts as a preservative to reduce or minimize growth of microorganisms within the dressing.

Mepilex Transfer Ag has been shown to inactivate microorganisms for up to 14 days *in vitro*.

Mepilex Transfer Ag consists of:

- a Safetac® soft silicone wound contact layer
- a compressed polyurethane foam containing silver sulfate colored by activated carbon

Safetac® is a proprietary adhesive technology developed by Mölnlycke that minimises pain to patients and trauma to intact skin or wounds.

Dressing material content:

Polyurethane, silicone, polyolefin, polyacrylate and silver sulfate which equals to 1.2 mg/ cm^2 silver

Indications for use

Mepilex Transfer Ag is indicated for the management of a wide range of exuding wounds such as leg and foot ulcers, pressure ulcers, partial thickness burns, traumatic and surgical wounds.

Mepilex Transfer Ag can also be used under compression bandaging.

Precautions

- Do not use Mepilex Transfer Ag on patient and/or user with known hypersensitivity to the ingoing materials/components of the product.
- Clinicians/healthcare professionals should be aware that there are very limited data on prolonged and repeated use of silver containing dressings, particularly in children and neonates.
- In the event of clinical infection Mepilex Transfer Ag does not replace the need for systemic therapy or other adequate infection treatment.
- Prior to commencing radiation therapy remove Mepilex Transfer Ag if product is present in the treatment area. A new dressing can be applied following treatment.
- Avoid contact with electrodes or conductive gels during electronic measurements, e.g. electrocardiograms (ECG) and electroencephalograms (EEG).
- Do not use Mepilex Transfer Ag together with oxidising agents such as hypochlorite solutions or hydrogen peroxide.
- Other than saline solution or water, the interaction of cleansing agents in combination with Mepilex Transfer Ag has not been demonstrated.
- The interaction of Mepilex Transfer Ag with topical treatments has not been demonstrated.
- Frequent or prolonged use of this product may result in transient or permanent discoloration of skin.
- For external use only.
- Do not reuse. If reused performance of the product may deteriorate, cross contamination may occur.
- Sterile. Do not use if sterile barrier is damaged or opened prior to use. Do not re-sterilise.
- Do not use after expiry date. If the product is used after the expiry date product properties cannot be ensured.

Instructions for use

Mepilex Transfer Ag should be used by or under the supervision of a qualified health care professionals.

Note that local hygiene procedures should be followed prior to and following the dressing change.

1. Cleanse the wound with saline solution or water according to standard clinical practice.
2. Dry the surrounding skin thoroughly.
3. Select appropriate dressing size. For best clinical outcomes, ensure that the product covers the wound bed and overlaps the wound edges. If more than one piece of Mepilex Transfer Ag is required, overlap the edges.
4. Cut to appropriate size (if needed). Remove the release films and apply the product with the adherent side to the wound. Do not stretch.
5. Apply a secondary dressing. Choose secondary dressing appropriate for the level of exudate. The size of the secondary dressing should exceed the edges of Mepilex Transfer Ag. Fixate to secure close contact with the wound.

Each Mepilex Transfer Ag dressing may be left in place for up to 14 days, depending on the patient, the condition of the wound and surrounding skin, or as indicated by accepted clinical practice.

Mepilex Transfer Ag is intended for short-term use up to 4 weeks. For long-term use, a clinical assessment by a physician is recommended.

Special storage conditions and handling conditions

Mepilex Transfer Ag should be stored in dry conditions below 25°C/77°F and protected from direct sunlight. Note that any color variances in Mepilex Transfer Ag do not affect product effectiveness or safety.

Disposal should be handled according to local environmental procedures.

Other information

If any serious incident has occurred in relation to the use of Mepilex Transfer Ag, you should report to Mölnlycke Health Care and to your local competent authority.

Mepilex and Safetac are registered trademarks of Mölnlycke Health Care AB.

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