Mölnlycke®

Mepilex® Border Flex Lite

WITH SAFETAC® TECHNOLOGY

FOAM DRESSING

Revised 2025-02

40681-45

Master PD-709689 Rev. 1

PD-753840 Rev. 2

Manufacturer / Fabricant
Mölnlycke Health Care AB
Entreprenörsstråket 21
SE-431 53 Mölndal, Sweden

NATURAL RUBBER LATEX


This product and its packaging are not made with natural rubber latex

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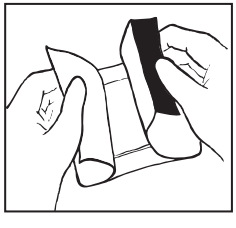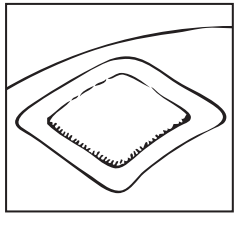
Made in Finland

What is Mepilex® Border Flex Lite? en

Mepilex Border Flex Lite is a flexible and conformable self-adherent, absorbent and showerproof dressing that maintains a moist wound environment. The waterproof backing layer protects the wound from dirt and bacteria. Safetac Innovative Adhesive Technology, minimizes pain when removing dressings, and minimizes damage to wounds and surrounding skin when removing dressings.

About Safetac® Technology Adhesive Layer

Safetac Technology Adhesive Layer keeps the dressing in place while gentle to skin and easy to use.

You can find Safetac® Technology in the following wound care brands: Mepilex® Lite, Mepilex® Border Flex, Mepilex® Border Flex Lite, Mepilex® Border Post-Op and Mepitac®. Read more at www.molnlycke.com.

What is Mepilex® Border Flex Lite for?

Mepilex Border Flex Lite is designed for a range of light draining wounds, such as abrasions, blisters and minor surgical incisions.

Material content

Dressing: Silicone, polyurethane, polyacrylate, viscose, polyester, polyolefin.

Precautions

- Do not use if you or the patient have a known hypersensitivity to the dressing or its ingoing materials, see Material content
- Do not use together with oxidizing cleaning solutions such as hypochlorite solutions or hydrogen peroxide
- If you see signs of infection e.g. fever or the wound or surrounding skin becoming red, warm or swollen, consult your health care professional, nurse or doctor
- Do not reuse. If reused performance of the product may deteriorate, cross contamination may occur
- Sterile dressing. Do not use if inner package is damaged or opened prior to use
- If the status of the wound has not significantly improved within 7 days, consult your health care professional, nurse or doctor

Instructions for use

1. Clean the wound with mild soap and water
2. Dry the surrounding skin thoroughly
3. To ensure proper adhesion, do not use moisturizers on the skin covered by the dressing
4. Select an appropriate size. The wound pad should cover the dry surrounding skin by at least ½ inch
5. Remove the first protective film covering the adhesive side of the dressing
6. Apply the adhesive side to the wound
7. Remove the remaining protective film and smooth down the border on the skin. Do not stretch

The dressing change interval may be several days. Change the dressing before it is fully saturated, at signs of leakage, or as indicated by healthcare provider.

Special storage- and handling conditions

Should be stored in dry conditions. The foam may become more yellow when exposed to light, air and/or heat. This has no influence on product properties. Disposal should be handled according to local environmental procedures.

Other information

Mepilex Border Flex Lite is a single use product.

This product and its packaging are not made with Natural Rubber Latex (NRL).

If the product is used after expiry date, product properties cannot be ensured.

The Mölnlycke trademarks, including Mölnlycke®, Mepilex® and Safetac® name and logo are registered globally to one or more of the Mölnlycke Health Care group of companies.

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