Mode d'emploi

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
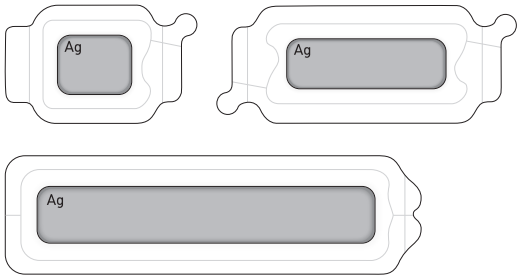
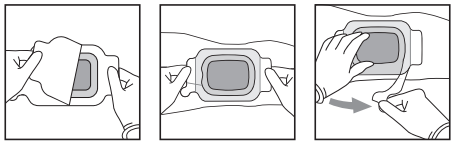
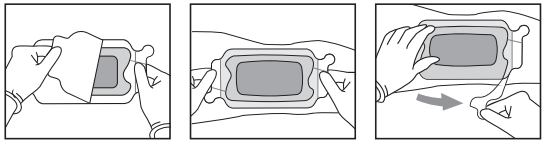
For product 10×35 cm

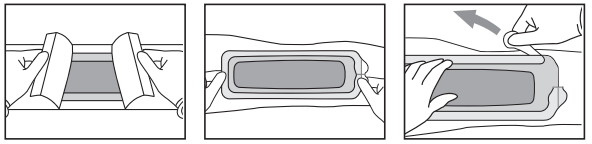en Mepilex® Border Post-Op Ag

Self-adherent antimicrobial soft silicone foam dressing

℞ only

Caution: Federal [USA] Law restricts this device to sale by or on the order of a physician [or properly licensed practitioner].

Product description

Mepilex Border Post-Op Ag consists of a Safetac® soft silicone wound contact layer, an absorbent polyurethane foam pad containing a silver compound and activated carbon, a layer with super absorbent polyacrylate fibres, a non woven and a vapour permeable waterproof film.

Mepilex Border Post-Op Ag is a soft silicone foam dressing that absorbs exudate, maintains a moist wound healing environment and has antimicrobial properties.

Mepilex Border Post-Op Ag has been shown to inactivate wound related pathogens up to 7 days in vitro.

Dressing material content:

Polyurethane, silicone, polyacrylate, cotton, viscose, polyester, polyolefin, silver sulfate which equals to 1.2 mg/cm^2 silver.

Indications for use

Mepilex Border Post-Op Ag dressing is indicated for the management of exuding wounds such as surgical and traumatic wounds.

Mepilex Border Post-Op Ag can also be used under compression bandaging.

Silver sulfate present in the dressing helps reduce microbial colonization in the dressing.

Precaution(s)

- Do not use Mepilex Border Post-Op Ag on patient and/or user with known hypersensitivity to the ingoing materials/components of the product.
- Clinicians/healthcare professionals should be aware that there are very limited data on prolonged and repeated use of silver containing dressings, particularly in children and neonates.
- Mepilex Border Post-Op Ag may cause transient discoloration of the wound bed and surrounding skin. Frequent or prolonged use of this product may result in permanent discoloration of skin.
- In the event of clinical infection Mepilex Border Post-Op Ag does not replace the need for systemic therapy or other adequate infection treatment.
- Prior to commencing radiation therapy remove Mepilex Border Post-Op Ag if product is present in the treatment area. A new dressing can be applied following treatment.
- Avoid contact with electrodes or conductive gels during electronic measurements, e.g. electrocardiograms (ECG) and electroencephalograms (EEG).
- Do not use Mepilex Border Post-Op Ag together with oxidising agents such as hypochlorite solutions or hydrogen peroxide.
- Other than saline solution or water, the interaction of cleansing agents in combination with Mepilex Border Post-Op Ag has not been demonstrated.
- The interaction of Mepilex Border Post-Op Ag with other topical treatments has not been demonstrated.
- For external use only.
- Do not reuse. If reused performance of the product may deteriorate, cross contamination may occur.
- Sterile. Do not use if sterile barrier is damaged or opened prior to use. Do not re-sterilise.
- Do not use after expiry date. If the product is used after the expiry date product properties cannot be ensured.

Instructions for use

Mepilex Border Post-Op Ag should be used by or under the supervision of a qualified health care professional.

Note that local hygiene procedures should be followed prior to and following the dressing change.

1. Cleanse the wound with saline solution or water according to standard clinical practice.
2. Dry the surrounding skin thoroughly.
3. The Mepilex Border Post-Op Ag wound pad should overlap the wound by at least 1-2 cm in order to protect the surrounding skin from maceration and fixate the dressing securely.
4. For dressings with one release film, remove the release film partly or completely. For dressings with two release films, remove the small release film partly or completely and then remove the large release film, partly or completely.
5. Apply the adherent side to the wound, allowing initial fixation of the dressing. If applicable, gently remove the remainder of the release film(s) while smoothening the dressing with the other hand to avoid wrinkles. Do not stretch the dressing during application.
6. Grip the application frame and tear the perforation apart.
7. Gently remove the application frame. Avoid wrinkles while fixating the edges by smoothening the dressing with the other hand.

Mepilex Border Post-Op Ag may be left in place for up to seven days, depending on the patient, the condition of the wound and surrounding skin, or as indicated by accepted clinical practice.

Initially, it might be required to change Mepilex Border Post-Op Ag more frequent due to that a change in treatment regime may result in an initial increase in exudation.

Mepilex Border Post-Op Ag is intended for short-term use up to 4 weeks. For long-term use, a clinical assessment by a physician is recommended.

Special storage conditions and handling conditions

Mepilex Border Post-Op Ag should be stored in dry conditions below 25 °C/77 °F and protected from direct sunlight. Note that any color variances in Mepilex Border Post-Op Ag do not affect product performance or safety.

Disposal should be handled according to local environmental procedures.

Other information

If any serious incident has occurred in relation to the use of Mepilex Border Post-Op Ag, it should be reported to Mölnlycke Health Care and to your local competent authority.

Mepilex® and Safetac® technology are registered trademarks of Mölnlycke Health Care AB.

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